# abcam

### Product datasheet

## Mouse AIF ELISA Kit ab235651

Recombinant SimpleStep ELISA

#### 1 References 3 Images

#### Overview

**Precision** 

Recovery

Mouse AIF ELISA Kit **Product name** 

**Detection method** Colorimetric

Sample	n	Mean	SD	CV%	
Cell extract	8			5%	

Sample type Cell culture extracts

Assay type Sandwich (quantitative)

Sensitivity 5.9 ng/ml

7.81 ng/ml - 500 ng/ml Range

Sample type	Average %	Range
Cell culture extracts	86	% - %

Assay time 1h 30m

**Assay duration** One step assay

**Species reactivity** Reacts with: Mouse

**Product overview** Mouse AIF ELISA Kit (ab235651) is a single-wash 90 min sandwich ELISA designed for the

quantitative measurement of AIF protein in cell culture extracts. It uses our proprietary SimpleStep

ELISA® technology. Quantitate Mouse AIF with 5.9 ng/ml sensitivity.

SimpleStep ELISA® technology employs capture antibodies conjugated to an affinity tag that is recognized by the monoclonal antibody used to coat our SimpleStep ELISA® plates. This approach to sandwich ELISA allows the formation of the antibody-analyte sandwich complex in a single step, significantly reducing assay time. See the SimpleStep ELISA® protocol summary in the image section for further details. Our SimpleStep ELISA® technology provides several benefits:

1

Intra-assay

Sample specific recovery

- Single-wash protocol reduces assay time to 90 minutes or less
- High sensitivity, specificity and reproducibility from superior antibodies
- Fully validated in biological samples
- 96-wells plate breakable into 12 x 8 wells strips

A 384-well SimpleStep ELISA® microplate (ab203359) is available to use as an alternative to the 96-well microplate provided with SimpleStep ELISA® kits.

**Notes** 

Abcam has not and does not intend to apply for the REACH Authorisation of customers' uses of products that contain European Authorisation list (Annex XIV) substances.

It is the responsibility of our customers to check the necessity of application of REACH Authorisation, and any other relevant authorisations, for their intended uses.

**Platform** 

Pre-coated microplate (12 x 8 well strips)

#### **Properties**

#### Storage instructions

Store at +4°C. Please refer to protocols.

Components	1 x 96 tests
10X Mouse AIF Capture Antibody	1 x 600µl
10X Mouse AIF Detector Antibody	1 x 600µl
10X Wash Buffer PT (ab206977)	1 x 20ml
50X Cell Extraction Enhancer Solution (ab193971)	1 x 1ml
5X Cell Extraction Buffer PTR (ab193970)	1 x 10ml
Antibody Diluent 4BR	1 x 6ml
Mouse AIF Lyophilized Recombinant Protein	2 vials
Plate Seals	1 unit
Sample Diluent NS (ab193972)	1 x 12ml
SimpleStep Pre-Coated 96-Well Microplate (ab206978)	1 unit
Stop Solution	1 x 12ml
TMB Development Solution	1 x 12ml

#### Function

Probable oxidoreductase that has a dual role in controlling cellular life and death; during apoptosis, it is translocated from the mitochondria to the nucleus to function as a proapoptotic factor in a caspase-independent pathway, while in normal mitochondria, it functions as an antiapoptotic factor via its oxidoreductase activity. The soluble form (AlFsol) found in the nucleus induces 'parthanatos' i.e., caspase-independent fragmentation of chromosomal DNA. Interacts with ElF3G,and thereby inhibits the ElF3 machinery and protein synthesis, and activates casapse-7 to amplify apoptosis. Plays a critical role in caspase-independent, pyknotic cell death in

#### Involvement in disease

hydrogen peroxide-exposed cells. Binds to DNA in a sequence-independent manner.

Defects in AIFM1 are the cause of combined oxidative phosphorylation deficiency type 6 (COXPD6) [MIM:300816]. It is a mitochondrial disease resulting in a neurodegenerative disorder characterized by psychomotor delay, hypotonia, areflexia, muscle weakness and wasting.

Sequence similarities

Post-translational

modifications


Belongs to the FAD-dependent oxidoreductase family.

Under normal conditions, a 54-residue N-terminal segment is first proteolytically removed during or just after translocation into the mitochondrial intermembrane space (IMS) by the mitochondrial processing peptidase (MPP) to form the inner-membrane-anchored mature form (AlFmit). During apoptosis, it is further proteolytically processed at amino-acid position 101 leading to the generation of the mature form, which is confined to the mitochondrial IMS in a soluble form (AlFsol). AlFsol is released to the cytoplasm in response to specific death signals, and translocated to the nucleus, where it induces nuclear apoptosis in a caspase-independent manner.


**Cellular localization** 

Mitochondrion intermembrane space. Mitochondrion inner membrane. Cytoplasm. Nucleus. Cytoplasm > perinuclear region. Proteolytic cleavage during or just after translocation into the mitochondrial intermembrane space (IMS) results in the formation of an inner-membrane-anchored mature form (AlFmit). During apoptosis, further proteolytic processing leads to a mature form, which is confined to the mitochondrial IMS in a soluble form (AlFsol). AlFsol is released to the cytoplasm in response to specific death signals, and translocated to the nucleus, where it induces nuclear apoptosis. Colocalizes with ElF3G in the nucleus and perinuclear region.

#### **Images**

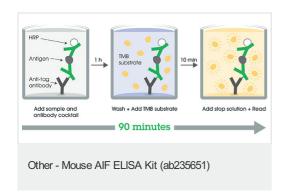


SimpleStep ELISA technology allows the formation of the antibodyantigen complex in one single step, reducing assay time to 90 minutes. Add samples or standards and antibody mix to wells all at once, incubate, wash, and add your final substrate. See protocol for a detailed step-by-step guide.

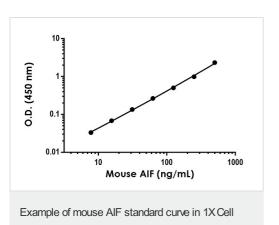

Extraction Buffer PTR.

Background-subtracted data values (mean +/- SD) are graphed.

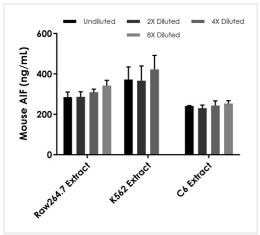

Interpolated concentrations of native AIF in Raw264.7 cell extract, K562 cell extract, and C6 cell extract based on a 75 µg/mL, 225 µg/mL, and 250 µg/mL extract load, respectively.

The concentrations of AIF were measured in duplicate and interpolated from the AIF standard curve and corrected for sample dilution. The interpolated dilution factor corrected values are plotted (mean +/- SD, n=2). The mean AIF concentration was determined to be 305.9 ng/mL in Raw264.7 cell extract, 387.4 ng/mL in K562 cell extract, and 242.1 ng/mL in C6 cell extract.

Please note: All products are "FOR RESEARCH USE ONLY. NOT FOR USE IN DIAGNOSTIC PROCEDURES"

#### Our Abpromise to you: Quality guaranteed and expert technical support

- · Replacement or refund for products not performing as stated on the datasheet
- · Valid for 12 months from date of delivery
- Response to your inquiry within 24 hours
- We provide support in Chinese, English, French, German, Japanese and Spanish
- Extensive multi-media technical resources to help you
- We investigate all quality concerns to ensure our products perform to the highest standards

If the product does not perform as described on this datasheet, we will offer a refund or replacement. For full details of the Abpromise, please visit https://www.abcam.com/abpromise or contact our technical team.

#### Terms and conditions

· Guarantee only valid for products bought direct from Abcam or one of our authorized distributors